

U2 Wrist Brace and Wrist Brace with Thumb

NOTE: Instructions cover multiple product families, images may vary.



Figure 1



Figure 2

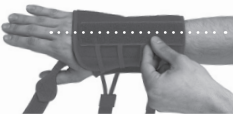


Figure 3



Figure 4



Figure 5



Figure 6

**NOT MADE WITH
natural rubber latex**

Please read the following instructions carefully and completely before use. Correct application is required to ensure proper function of the device.

Intended for single patient use.



WARNING: If your pain increases or persists or you experience any unusual reactions while using this product, consult your healthcare provider immediately.

APPLICATION INSTRUCTIONS

1. Unfasten straps and slide affected hand into brace (Figure 1).
2. Position palm stay firmly against the palm of the hand (Figure 2). Ensure the edge of the brace is below the knuckles (Figure 3).
3. Detach and align dorsal stay with middle knuckle. Reattach (Figure 3).
4. Secure the thumb-web strap (Figure 4).
5. Pull and wrap the lacing straps around wrist. Secure when the desired fit and level of compression is reached (Figure 6).
6. Adjust straps as needed when swelling of the hand, wrist or thumb changes.

NOTE: Pediatric braces do not have a detachable dorsal stay

NOTE: For braces with thumb support; wrap thumb strap around thumb and secure (Figure 5).

NOTE: Please refer to sew-in label on device for care instructions.

Please visit our website at www.birdcronin.com for additional product information or contact a customer care representative at **800-328-1095**.

NOTE: Every effort has been made to obtain the maximum compatibility of function, strength, durability and comfort with this product, however, there is no guarantee that injury will be prevented with its use.